

CBX portable balances are ideal for use in education, laboratories, offices or production control. The bright LCD display shows 18mm digits and includes a low battery indicator. No need to replace batteries as the CBX comes standard with an internal rechargeable battery. The battery can last up to 216 hours before needing a recharge, which can be done via a low-voltage USB-type input. The CBX can also operate with a power adapter, which is included with the balance. The easy-to-clean 304 stainless steel tray protects the CBX from most materials and liquids. Calibrating the CBX is easily done with the keypad to ensure accurate and precise results at all times.



CBX balances include a standard rechargeable battery for a complete portable weighing solution



Features

- Solid ABS housing is durable for industrial use
- External calibration allows for verification and adjustment with weights
- Vivid, backlit LCD easily visible in any lighting conditions
- Auto power-off to save energy
- Simple key operation allows for easy set up and operation
- Large, grade 304 stainless steel pan allows swift cleaning
- Level indicator and adjustable feet ensure proper balance setup for optimum weighing results
- Rechargeable battery included for operation almost anywhere
- Low-battery indicator alerts users when batteries need replacement
- Parts counting with preset sample sizes
- Operates with a low voltage USB input providing flexibility for powering the device

Applications


- Weighing
- Parts counting

About ADAM:

For nearly 50 years, professionals worldwide have relied on Adam Equipment for an extensive selection of dependable and affordable weighing equipment. Adam provides the right balance of speed, performance and value.

Product Details

Model	CBX 1201	CBX 3001	CBX 6001
Capacity	1200g	3000g	6000g
Readability	0.1g	0.5g	1g
Repeatability	0.2g	1g	2g
Linearity	0.3g	1.5g	3g
Pan Size	5.8x5.8" / 148x148mm		
Weighing Units	g, N, lb, oz		
Calibration	External		
Display	Backlit LCD with 0.7" / 18mm high digits		
Power Supply	Internal rechargeable battery / Power Cord		
Operating Temperature	0° to 40°C / 32° to 104°F		
Housing	ABS plastic		
Overall Dimensions	6.1"x8.9"x2.4" / 155x225x60mm		
Net Weight	1.3lb / 0.572kg		

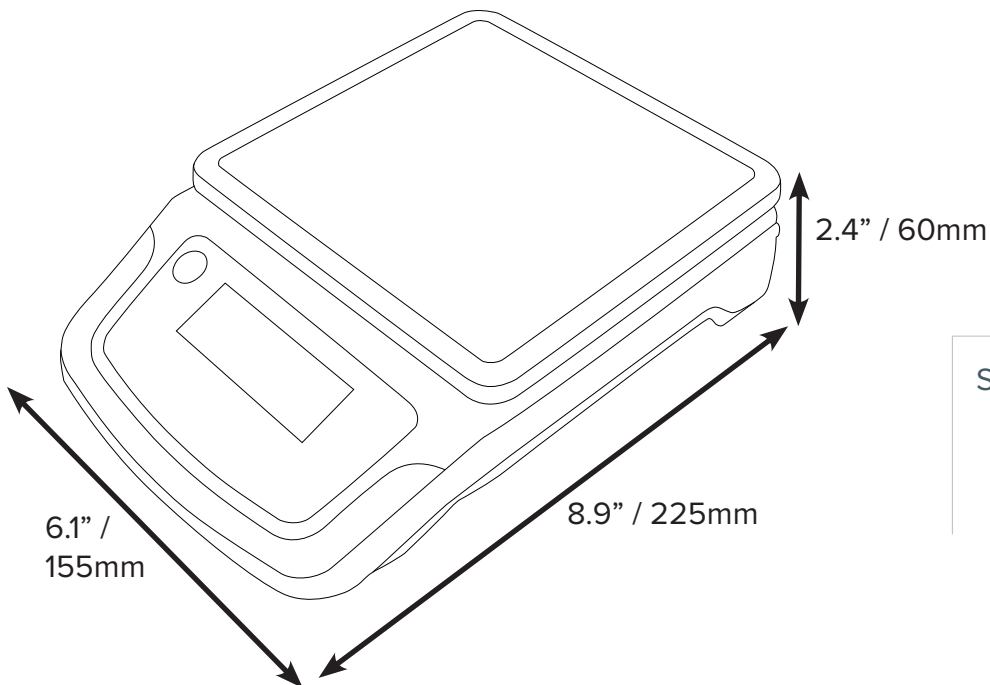

Supplied by:

www.adamequipment.com

Adam Equipment follows a policy of continuous improvement and reserves the right to change specification and pricing without notice. Adam Equipment's standard terms and conditions of sale apply. Errors and omissions excepted.